1 Warning

- Do not use antiperspirant immediately after irradiating the armpit.
- Please confirm and clean the dirt on the irradiation port regularly to avoid burns and malfunctions during use.
- Do not use it with the vent blocked, to avoid smoke and fire.
- Do not point the irradiation port towards a black object, so as to avoid the machine not working normally.
- Do not irradiate clothes, hair and other items to avoid charring, discoloration and damage.
- Do not use this product around the eyes to avoid eye damage.
- After use, please wipe off the hair and skin care products attached to the epilator to avoid malfunction and skin damage.
- Do not use around children, the hair removal device will emit strong light, which can easily hurt the eyes.
- Please set the irradiation intensity according to your own skin condition, and please start from the gentle mode for the first time
- It is normal for the skin to become slightly hot and red after use.
- The product cannot be left unattended after it is turned on, please turn it off in time after use.

- Clean the skin to ensure that the skin surface is free of hair residue, completely dry and free of oily substances and stains.
- Do not expose the product to sunlight or UV light for a long time, which will seriously damage the product.
- If the built-in fan does not run after the power is turned on (no humming sound is heard), do not use it to avoid malfunctions such as smoke and fire.

Unsuitable people and situations

- Do not use for pregnant women, breastfeeding mother with a child under six months old, women's menstrual period, body surface scars, diabetes, fragile skin and light sensitive skin.
- Severe heart disease, moderate to severe hypertension, malignant tumor, infectious skin disease, burnt skin and anticoagulant disorders are prohibited.
- Do not use if you have abnormal collagen metabolism, including a history of scarring and poor wound healing.
- People with vascular disease, such as varicose veins and vasodilation at the irradiation site, do not use it.
- Do not use if the skin is sensitive to light and prone to rashes or allergies.
- It is contraindicated in any hemorrhagic patient, and external use of isotretinoin, fruit acid (AHA), salicylic acid (BHA) or azelaic acid is contraindicated.
- Those who use chemical hair removal within 30 days, those who have used skin whitening products in a short period of time, or those who have used cosmetics containing vitamin A, please use it with caution.
- It is forbidden to take aspirin drugs, immunosuppressants, isotretinoin and other drugs, and those who have recently used light-sensitive drugs (such as tetracycline and sulfonamide, etc.).
- Can not be used at the same time as essential oil products, please ensure that the skin of the irradiated area is clean and dry, and no lotions, lubricants and other cosmetics are used.

●Not applicable body parts

1/ Eyeballs, eyelids, around the eyes 2/ Adam's apple 3/ Inside of mouth, genitals, mucous membranes, etc. 4/ Plastic surgery, surgical site 5/ Nipple, areola 6/ Spots, moles 7/ Eyebrows, hair 8/ Inside of navel 9/ Acne Acne area 10/ Tattoo, body painting area 11/Head 12/Large vein 13/Ear (\*Note: This list is not exhaustive, if you are not sure whether you can use this product, it is recommended that you consult a doctor)

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Warning: Do not use this product near water-filled bathtubs, shower trays, washbasins or other water-filled utensils, keep this product and power adapter dry.

- If the product is damaged, do not touch any internal parts to avoid electric shock. Do not insert any object into this product, and do not place this product at will.
- Do not use this product around the eyes (even if you wear protective glasses), if you operate around the eyes, the light will be absorbed by the iris melanin in your eyes at close range, resulting in danger.
- In order to avoid the impact of intense pulsed light on your vision, please wear protective glasses during each operation.
- It is recommended that you seek help from others when doing back hair removal to avoid danger
- In order to prevent the excessively long hair from absorbing heat on the skin surface, causing redness and burning of the skin and odor, if your hair is longer than 1mm, please use a scraper to shave it before operating the machine.
- In order to ensure the use effect and avoid the product from not flashing or accidentally flashing, please ensure that the depilation area is dry and clean.
- Due to the lack of data on the use of this product on children under the age of 15, please pay attention to preventing children under the age of 15 from contacting this product, and please keep away from children under the age of 15 when operating the instrument.
- To ensure the effect, please use different light intensity modes for different body parts as needed.
- If the power adapter is damaged, in order to avoid danger, it must be replaced by professionals in the manufacturer's maintenance department or similar departments.

▲Precautions

- 1-Do not use this instrument for naturally dark skin;
- 2-Do not expose to the sun before hair removal;
- 3- It is forbidden to use it in a dark environment;
- 4-Do not hit and shake the product violently;
- 5-Do not use in direct sunlight to avoid machine failure;
- 6- Please wear protective glasses when using;
- 7-Please shave before hair removal (do not use chemical hair removal method);
- 8- Be sure to make the hair removal window close to the skin when removing hair, and then perform the hair removal operation;
- 9-Be sure to start from the gentle mode when using it for the first time;
- 10-Do not pull out the plug directly when the machine is not turned off to avoid damage to the machine.

About Sapphire Ice Cooling Hair Removal Device

·This product relies on the concave thermal flat hair removal device to miniaturize the medical beauty hair removal technology and make it home-use. During hair removal, the flat head compresses the skin and overwhelms the hair follicles, making the roots closer to the higher energy area of the epidermis and improving the energy absorption efficiency of the hair follicles.


·Advanced Sapphire ice-cooling system maintains 5°C ice contact during the treatment. More chilling, more comfortable.

·The flat lamp head is easy to clean, clean after wiping, with stronger light transmission and better hair removal effect.

·This product has 3 gears adjustment to meet the needs of different skin sensitivities and is more humane. The fuselage and the power supply are detachable, compact and portable, and easy to store. This product has a high-life, one-piece lamp head that does not need to be replaced, and is suitable for multiple hair removal areas (such as: face, lips, underarms, bikini line, arms, legs, etc.) to save you expensive lamp replacement costs, No frequent accessory purchases and no hidden additional fees. In addition, it has a perfect filtering system, which can well filter invalid short wavelength and harmful light.

This product comes with LED red light skin rejuvenation mode. After it is turned on, the instrument can continuously use red light to irradiate the skin, which has a certain effect on the skin.

Product structure

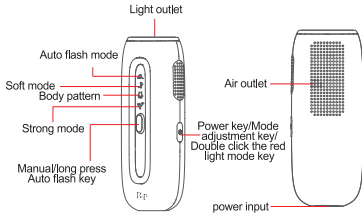●Preparation before use

Clean skin thoroughly. Make sure the skin is hair-free, completely dry, and free of oily residue. (residues such as perfumes, lotions, body lotions and sunscreens)Make sure that there is no dirt, hair, etc. remaining in the irradiation port, and that there is no damage or abnormality. If there is dirt, be sure to clean it up to avoid product failure, deterioration and skin problems. Perform a skin test before using the device, and confirm the skin condition 2 hours after the skin test, if there is no abnormality (redness, blisters) on the irradiated area, you can use this product with confidence.

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SKIN TONE						
HAIR COLOUR						
White/Grey	✗	✗	✗	✗	✗	✗
Light Blonde	✗	✗	✗	✗	✗	✗
Red	✗	✗	✗	✗	✗	✗
Dark Blonde	✓	✓	✓	✓	✗	✗
Light Brown	✓	✓	✓	✓	✗	✗
Brown	✓	✓	✓	✓	✗	✗
Dark Brown	✓	✓	✓	✓	✗	✗
Black	✓	✓	✓	✓	✗	✗
✓ It's match!	✗ Try a different hair removal method					

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## Instructions


●Step 1

▲ Do not touch the squeegee directly with the skin, use it with foam or shave in the shower. The scraper is sharp, please use it carefully to avoid scratches.

●Step 2

Note: After finishing the hair removal operation, be sure to turn off the power before unplugging the power cord to prevent damage to the power supply.

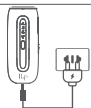---

● Step 3

**Note:** Do not look directly at the flash when removing hair.

●Step 4 

Long press the power button to turn on, short press the power button to adjust the mode according to your needs, double-click the power button after booting to turn on the red light skin rejuvenation mode, and double-click again to turn off the red light skin rejuvenation mode.

●Step 5

According to your own needs, short press the manual flash button to emit light, the lamp head does not touch the skin, long press the flash button for 2S, the automatic flash mode indicator will light up, indicating that the automatic flash mode is turned on.

●Step 6

Place the light outlet on the area to be treated on the skin and fit it vertically. The indicator icon on the screen mode of the hair meter flashes, indicating that it has been fully contacted. Press the flash button to start the hair removal operation. If the indicator icon does not flash, it indicates that the lamp head and the skin are not fully fitted, and the instrument will not emit light.

▲ After the hair removal is complete, press and hold the power button to turn off the device (do not unplug the device without turning it off to avoid damage)

(suitable for large areas such as arms and legs)

Automatic flash mode: Long Press the flash button for 2s, the automatic flash mode indicator icon will light up, place the light output port on the skin to be treated area. The flash will be emitted at certain intervals automatically. Press and hold the flash button for 2s to turn off the automatic mode.

(suitable for small areas such as fingers and lip hair)

**Manual Flash Mode:** Place the light outlet in the skin to be treated area and fully fit the skin. Short press the flash button to start hair removal.

\*TIP: At any point during the hair removal process, if your skin starts to feel overheated, please pause the operation and allow the light outlet to cool down again before hair removal.

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Use 2 times a week.

Note: Please use it multiple times according to the growth of the hair. If the hair is too thick, the gear and frequency can be appropriately increased on this basis. The specific situation varies from person to person.

Face - 3mins, Armpit - 5mins, Arm - 5mins, Bikini line - 5mins, Leg - 10mins.

## After-use care

### ●Precautions after hair removal

- It is recommended to avoid touching the irradiated area with water within 6 hours after the hair removal operation, and then moisturizing the irradiated area with moisturizing water.
- Do not use perfumes, masks and cosmetics containing scrub beads after skin depilation.
- Within 24 hours after the hair removal operation, please avoid UV exposure and apply sunscreen if necessary.
- Artificial tanning is not permitted during hair removal operations;
- If you need plastic surgery after hair removal (only the hair removal area), please consult your doctor.
- Consult a surgeon if you feel unwell.

### ●Cleaning and maintenance

.In order to ensure the best results and longest lifespan of the product, please clean the product before and after each hair removal to avoid long-term accumulation of hair and dust, if you do not clean the product properly will not have the desired effect.

- Before starting cleaning, the equipment must be turned off and unplugged to allow the machine to cool.
- Wipe and clean the light outlet with a soft, clean and dry cloth.
- Store the equipment in a dry and dust-free environment at -10-60°C
- Do not wash this product and accessories in the faucet and dishwasher.
- Do not use steel wool, abrasive cleaners (toothpaste, etc.) or aggressive liquids (gasoline, motor oil, acetone, etc.) to clean the product.

### ●Environmental protection

When disposing of the product, please do not throw it with the general household waste, but tell it to the designated recycling center, which is good for the environment.

## Basic parameters

Product Name	: Sapphire Ice Cooling IPL Hair Removal Device
Product Model	: W-1098
Product size	: 163*60*37mm
Rated input	: 12V 3A
Power	: 36W
Energy level	: 3 gears
Intensity	: 12-16J/ cm <sup>2</sup>
Lamp head area	: 3 cm <sup>2</sup>
Spectrum	: 530-1200nm
Executive standard	: GB4706.1-2005

Exclusively Distributed by R&P International Company Limited

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